



ALIBI Reader Software

ALIBI Memory data readout

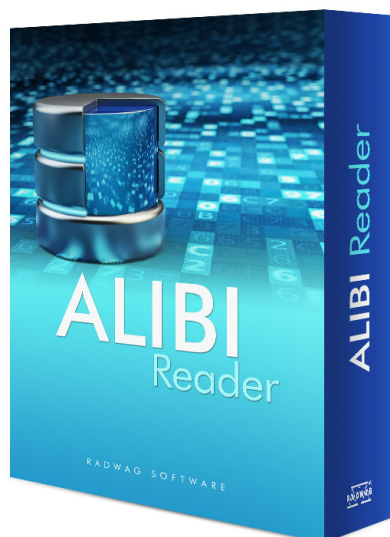

ALIBI Reader

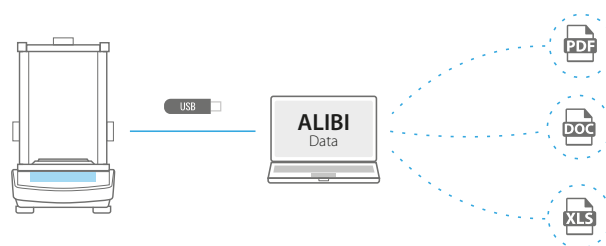

ALIBI Memory export

Features

Transfer and Configuration of Data

ALIBI Reader enables previewing previously exported file containing ALIBI memory data on a PC computer. The software allows to configure and filter data, e.g. by date, operator code, product, and saving the data as CSV file.


Generating Reports

ALIBI Reader facilitates creating reports of previously selected measurements and saving them as XLS, DOC and PDF files.

Compatibility with Weighing Instruments

The software is compatible with X2 and R series balances and X2 and R series moisture analyzers using PUE C32, PUE HX5.EX and PUE HX7 terminals. ALIBI Reader features five language versions: Polish, English, German, French and Spanish.

Screenshots

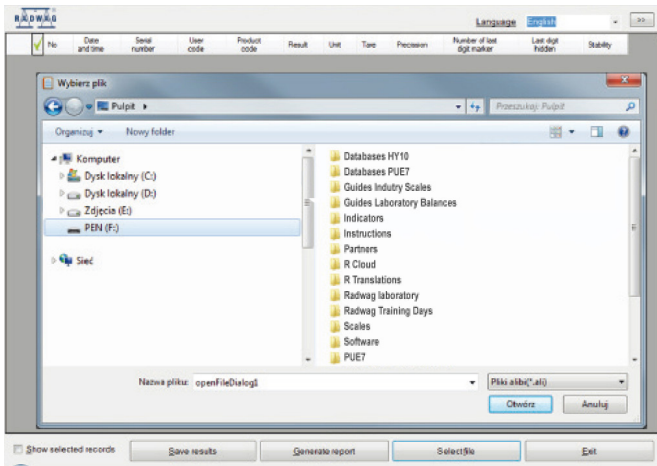

Import of the ALIBI file to ALIBI Reader

No	Date and time	Serial number	User code	Product code	Result	Unit	Tare	Precision	Number of last digit marker	Last digit hidden	Stability
1111	2014.01.03.12.19.10	12310000		U4624	13.1048	kgf1	0.0000	4	0	No	No
1112	2014.01.03.12.19.10	12310000		U4624	17.1261	sole	0.0000	4	0	No	No
1113	2014.01.03.12.19.10	12310000		U4624	17.1273	sole	0.0000	4	0	No	No
1114	2014.01.03.12.19.10	12310000		U4624	17.1273	sole	0.0000	4	0	No	Yes
1116	2014.01.03.12.19.20	12310000		U4624	17.1273	sole	0.0000	4	0	No	Yes
1116	2014.01.03.12.19.20	12310000		U4624	9.9885	pat1	0.0000	4	0	No	Yes
1117	2014.01.03.12.19.20	12310000		U4624	9.9885	pat1	0.0000	4	0	No	Yes
1112	2014.01.03.12.19.10	12310000		U4624	17.1261	sole	0.0000	4	0	No	No
1113	2014.01.03.12.19.10	12310000		U4624	17.1273	sole	0.0000	4	0	No	No
1114	2014.01.03.12.19.10	12310000		U4624	17.1273	sole	0.0000	4	0	No	Yes
1115	2014.01.03.12.19.20	12310000		U4624	17.1273	sole	0.0000	4	0	No	Yes
1116	2014.01.03.12.19.20	12310000		U4624	9.9885	pat1	0.0000	4	0	No	Yes
1117	2014.01.03.12.19.20	12310000		U4624	9.9885	pat1	0.0000	4	0	No	Yes
1128	2014.01.03.12.21.05	12310000	78GE	U4624	199.79	g	0.00	2	0	No	No
1129	2014.01.03.12.21.05	12310000	78GE	U4624	199.79	g	0.00	2	0	No	No
1130	2014.01.03.12.21.05	12310000	78GE	U4624	199.79	g	0.00	2	0	No	No
1131	2014.01.03.12.21.10	12310000	78GE	U4624	199.74	g	0.00	2	0	No	No
1132	2014.01.03.12.21.12	12310000	78GE	U4624	0.19676	kg	0.00000	5	0	No	Yes
1133	2014.01.03.12.21.12	12310000	78GE	U4624	0.19676	kg	0.00000	5	0	No	Yes
1134	2014.01.03.12.21.13	12310000	78GE	U4624	0.19679	kg	0.00000	5	0	No	No
1135	2014.01.03.12.21.13	12310000	78GE	U4624	0.19679	kg	0.00000	5	0	No	No
1136	2014.01.03.12.21.13	12310000	78GE	U4624	0.19679	kg	0.00000	5	0	No	Yes

ALIBI memory weighings

Compatible with

Weighing Instruments

- X2 laboratory balances
- R laboratory balances
- X2 Moisture Analyzers
- R Moisture Analyzers

Indicators and Terminals

- PUE C32
- PUE HX5.EX
- PUE HX7

Hardware recommendations

PC computer

- operating system: Windows 10 / 8 / 8.1 / 7 (excluding Starter version)
- 2 GHz dual-core processor or faster
- minimum 2 GB operating memory
- minimum 10 GB of free space on a hard drive
- 1280x1024 computer monitor
- DVD drive
- default system font size: 96 dpi