

R Panel


The R Panel software enables operation of selected weighing instruments and electromagnetic weighing modules via a computer operation panel



R Panel



Weighing result window in R Panel software

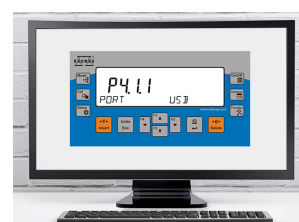

USB port selection screen

Features

Features

- Remote weighing instruments operation
- Creating screenshots
- Communication via RS232, USB and Ethernet
- Various languages
- Compatibility with Windows Vista, 7, 8, 8.1, 10, Server 2008R2, 2012, 2016

Intended Use

The R Panel software enables operation of R and WTC series balance, C315 indicator-based weighing instruments and MAS R and MPS R weighing modules via a computer operation panel. The software is intended for computers based on Windows operating system. The operator has access to all keys and functions that are to be found on an operation panel, such as zeroing, taring, setup etc.

Auto Searching for Devices

An innovative algorithm for automatic search of weighing instruments and electromagnetic weighing modules is a guarantee for a remarkable ease of connection of a device to a computer.

Accessibility in R-Lab

Functionality of computer operation panel has been transferred to RLab PC software, this enables operation of weighing instruments and electromagnetic weighing modules without the necessity to run any additional software.

Compatible with

Weighing Instruments

- R laboratory balances
- WTC laboratory balances
- C315 industrial scales

Weighing modules

- MAS R weighing module
- MPS R weighing module

Hardware requirements

PC Computer

- Operating system: Windows 10 / 8 / 8.1 / 7 (excluding Starter version)
- 2 GHz dual-core processor or faster
- Minimum 3 GB of operating memory
- Minimum 35 GB of free space on a hard drive
- 1280x1024 computer monitor
- Default system font size: 96 dpi